Data Sheet

Manual Force Test Stands ES10 & ES20

32-1047 REV 1020

Page 1 of 2

These compact and economical test stands are ideal for low force applications with relatively small sample sizes. The ES10 and ES20 produce both tensile and compressive force up to 100 IbF (500 N). Choose between lever operation for quick tests or hand wheel operation for better control of motion. Combine these test stands with Mark-10 digital force gauges and gripping fixtures to create complete testing solutions.

Specifications

Load capacity:	100 lbF [500 N]
Maximum travel:	9" [229 mm]
Loading method:	
ES10:	Lever
ES20:	Hand wheel
Travel rate:	
ES10:	1.05 in [26.7 mm] per lever rev.
ES20:	0.083 in [2.1 mm] per hand wheel rev.
Weight:	10 lb [4.5 kg]
Warranty:	3 years (see individual statement for details)

Optional Equipment

Dial travel indicator ES001	1" displacement, 0.001" resolution. Mounts onto the test stand and can be repositioned.
Digital travel indicator ES002	$1^{\prime\prime}$ / 25 mm displacement, 0.0005" / 0.01 mm resolution. Includes zero and in/mm selector buttons. Mounts onto the test stand and can be repositioned.
Indicator mounting kits ES001NI (-8MM)	For customer-supplied indicators with 8 mm or 9.5 mm dia. stem sizes.
Base plate AC1060	Steel base plate has matrix of #10-32 threaded holes for fixture mount- ing





Dial travel indicator

Digital travel indicator





Ordering Information


Model No.	Description
ES10	Test stand, lever-operated, 100 lbF
ES20	Test stand, hand wheel-operated, 100 lbF
ES001	Dial indicator kit for ES10 / ES20, 1" travel
ES002	Digital indicator kit for ES10 / ES20, 1" / 25 mm travel
ES001NI	Mounting kit for 9.5 mm dia. indicator stem (indicator not included)
ES001NI-8MM	Mounting kit for 8 mm dia. indicator stem (indicator not included)
AC1060	Base plate, threaded holes



The ES20 is shown in a peel testing application with Series 5 digital force gauge, and film & paper grips.

In The Box



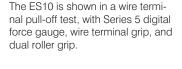



The ES10 and ES20 test stands are shipped as shown above, and include the following accessories:

- (1) Lever (ES10) or hand wheel (ES20)
- (1) Allen wrench set
- (4) Thumb screw for force gauge



Mark-10 Corporation www.mark-10.com info@mark-10.com Toll-free: 888-MARK-TEN I Tel: 631-842-9200 Fax: 631-842-9201



Specifications subject to change without prior notice


Data Sheet

Manual Force Test Stands ES10 & ES20

32-1047 REV 1020

Dimensions

in [mm]

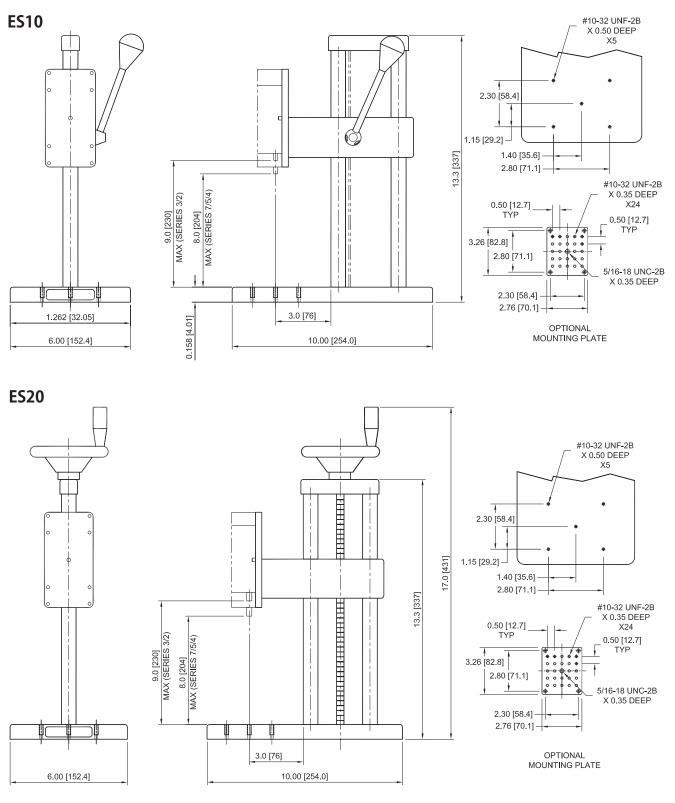



Mark-10 Corporation ■ www.mark-10.com ■ info@mark-10.com Toll-free: 888-MARK-TEN ■ Tel: 631-842-9200 ■ Fax: 631-842-9201



Page 2 of 2

Specifications subject to change without prior notice